


Transform STM PT

Compatible with ProTaper Gold

Shape with S files. Finish to the measured apical size.
ENDOTECH ROTARY · 400-500 rpm · 2 N·cm · MicroPath 13/.03: 1 N·cm

ENTRY CHECK Establish straight-line access and a reproducible hand-file path before rotary instruments enter that region.

THE INSTRUMENTS

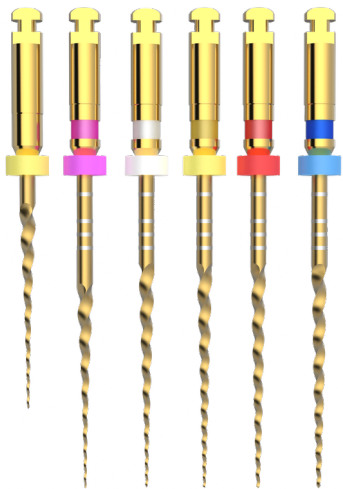

Original EndoTech instruments, left to right: SX, S1, S2, F1, F2, F3. Image is not to scale.

SEQUENCE + CLINICAL DECISIONS

- 1 Build the glide path**
Hand scout + MicroPath
Scout the coronal two-thirds with #10/#15 hand files. MicroPath 15/.03 or 17/.03 is the EndoTech glide-path option; add 13/.03 for complex or narrow paths.
- 2 Relieve coronal load**
S1 → S2; SX as needed
Use S1 then S2 only to the depth already negotiated. Brush on withdrawal, away from the furcation. Add SX for coronal relocation or flare.
- 3 Verify apical access**
#10 / #15 + WL
Negotiate the apical third. Establish WL, verify patency and confirm a loose, reproducible glide path. Extend the glide sequence if required.
- 4 Shape to length**
S1 → S2
Carry S1 then S2 to WL with controlled outstroke brushing. Reconfirm WL, irrigate and recapitulate before finishing.
- 5 Finish and gauge**
F1 → #20 hand gauge
Use F1 without brushing. Withdraw on reaching WL; do not dwell for more than one second. If the #20 hand file is snug at WL, assess the endpoint before enlarging.
- 6 Enlarge only if needed**
F2 / #25; F3 / #30
If the gauge is loose and enlargement is indicated, use F2 then gauge #25; use F3 then gauge #30 only if needed. Finish disinfection and cone fit.

GLIDE → S1 / S2 (+ SX) → F1 → GAUGE

SX: coronal work

SX is an auxiliary shaping instrument. Its place depends on access and coronal interference; it is not a compulsory full-length apical file. Preserve furcal dentin.

Two different movements

S1 / S2 / SX: controlled brushing on withdrawal. F1 / F2 / F3: passive, non-brushing progression with prompt withdrawal at WL. Do not dwell or repeatedly work the terminus.

Finishing is conditional

F2 and F3 are selected by gauging, anatomy and the preparation objective. Do not run the entire coloured row automatically. This card covers the six-file SX / S1 / S2 / F1 / F2 / F3 core.

Transform STM PT

Compatible with ProTaper Gold

Keep shaping and finishing distinct.

01 Familiar ProTaper sequence

This card adapts the ProTaper Gold shaping/finishing logic to the six EndoTech PT instruments. It preserves staged coronal preparation, apical glide verification, S1/S2 to length, then F1 and gauging-led enlargement.

03 Gauging is an assessment

Assess binding at WL, not resistance caused by coronal interference. Consider apical anatomy, curvature and the intended irrigation/obturation approach before selecting another finishing file.

02 SX placement


The original technique uses SX selectively for coronal relocation, shorter roots or additional flare. A visual SX-first product lineup is not an instruction to take SX to WL before the glide path.

04 Range and identification

S1 is purple, S2 white, F1 yellow, F2 red and F3 blue in the original artwork. F1/F2/F3 correspond to nominal 20/.07v, 25/.08v and 30/.09v. Variable taper labels do not describe a constant taper along the entire blade.

MOTOR SETUP

ROTARY: 400-500 rpm. PT SX, S1, S2, F1-F3 and MicroPath 15/.03 / 17/.03: 2 N-cm. MicroPath 13/.03: 1 N-cm.

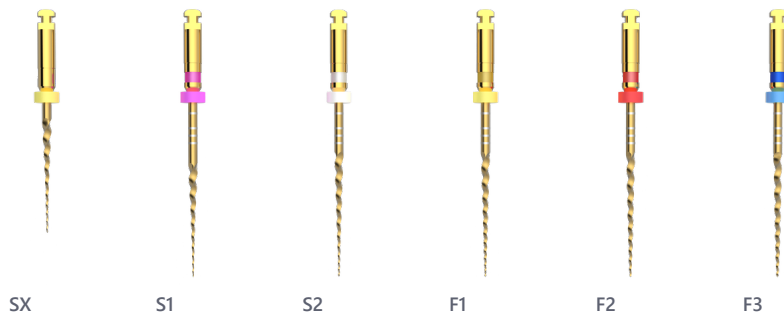

SOURCE TRAIL

[Dentsply Sirona | ProTaper Gold directions for use](#)

February 2017, pp. 2-3. Sequence, brushing distinction and apical gauging.

[Dentsply Sirona | ProTaper Gold 2025 IFU archive](#)

June 2025, English pp. 3-6. File designations, staged shaping, non-brushing finishing and gauging.

[EndoTech | Transform S PT product context](#)

Six-file source artwork and current local catalog. MicroPath steps are the EndoTech adaptation.