

Transform S™ RG

Compatible with WaveOne Gold

Start with Primary. Select another size from the canal assessment.

ENDOTECH MOTOR · Opener: 400-500 rpm / 2 N-cm · RG Glider + shaping: 170° CCW / 40° CW

ENTRY CHECK

Scout with a #08 or #10 hand file, establish WL and expand a reproducible glide path before full-length shaping.

THE INSTRUMENTS



Original EndoTech lineup: opener, Glider, Small, Primary, Medium, Large. Display order is not the clinical sequence.

SEQUENCE + CLINICAL DECISIONS

- 1 Prepare coronal entry**
17/.08 · 17 mm | ROTARY
Use the opener only in the known coronal path. Run its rotary program and brush away from the furcation. Stop before curvature and resistance.
- 2 Expand the glide path**
RG Glider | RECIPROCATING
After a loose #10 hand file reaches WL, use the dedicated Glider to WL. Confirm a reproducible path of at least ISO 15 and recheck WL.
- 3 Start shaping**
Primary 25/.07v
Begin with Primary in irrigant, using short, passive pecks. After 2-3 mm of advancement, withdraw, clean, irrigate, recapitulate and re-irrigate.
- 4 If Primary stalls**
Small 20/.07v → Primary
Reassess access and the glide path. If Primary still cannot progress, take Small passively to WL, then return to Primary. Never increase pressure.
- 5 If Primary is loose**
Medium 35/.06v; Large 45/.05v
If Primary is loose at WL with no apical dentin debris, assess whether Medium and then Large are indicated. Avoid automatic use of all four sizes.
- 6 Confirm and finish**
Endpoint → clean → seal
Remove promptly at WL. Assess apical engagement, anatomy and preparation adequacy; complete disinfection and verify cone fit under length control.

OPEN → GLIDE → PRIMARY → ASSESS

Primary is the start

Small is a branch for a Primary that cannot progress after the path is reassessed. The familiar WaveOne Gold logic begins shaping with Primary.

Read the response

Apical flute loading and looseness help guide the reference sequence. Debris alone does not prove adequate disinfection or an appropriate final size.

Switch the motion

Opener: continuous rotation. RG Glider and shaping: 170° counter-clockwise / 40° clockwise. If MicroPath is used instead of the RG Glider, switch to its rotary program.

Transform S™ RG

Compatible with WaveOne Gold

When to use Small, Medium or Large.

01 Reference logic

The WaveOne Gold guide starts with Primary, moves to Small if Primary will not advance, then returns to Primary. Medium or Large are selected when a larger apical preparation is justified. The opening step is EndoTech-specific.

03 Set the exact RG motion

Use 170° counter-clockwise / 40° clockwise for RG Glider and shaping files. The opener uses continuous rotation at 400-500 rpm and 2 N·cm. Check that the motor can deliver the exact angle pair before use.

02 Depth remains controlled


The same named file can progress in several passes. Coronal pre-enlargement is restricted to a region with a verified glide path. Reconfirm length after shaping changes the access path.

04 Identify the file

The lineup shows the opener, Glider, Small, Primary, Medium and Large. Check the designation on the pack before use. The Blender wireframe shows RG Primary.

MOTOR SETUP

RG GLIDER + SHAPING: 170° CCW / 40° CW. OPENER: rotary, 400-500 rpm; torque: 2 N·cm. CCW = counter-clockwise; CW = clockwise. Use a motor compatible with the specified motion.

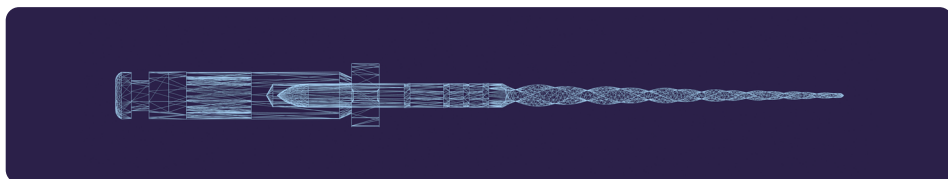

RG Primary · Native Blender wireframe render

SOURCE TRAIL

[Dentsply Sirona | WaveOne Gold directions for use](#)

March 2017, pp. 3-4. Glider entry, Primary-first logic and conditional size changes.

[EndoTech | Transform S RG product context](#)

Original opener / Glider / four-file artwork extracted from the EndoTech PT-RG brochure, p. 3.

[Dentsply Sirona | WaveOne Gold product resources](#)

Current manufacturer resource index; comparator context only.